

FS-1 series

Floor Stand



A versatile column and a stable baseplate



Topcon's CV-5000 refraction system compatible

The FS-1 floor stand consists of a versatile column and a stable baseplate, providing an economical solution suitable for the examination room as an alternative refraction lane. The FS-1 has a small footprint and accepts several models of phoropter holders such as a balanced phoropter arm or a parallel phoropter arm, as well as a chart projector bracket and a convenient reading light. The up-light fits seamlessly in the FS-1 design providing a homogeneous room light. The FS-1 can be combined with Topcon's OC-6 ophthalmic chair or one of the IS-1 series ophthalmic chairs, such as the OC-10 or OC-12. The baseplate with integrated wheels permits an easy displacement of the ophthalmic chair for easy wheelchair accessibility. The optional dark grey cabinet accommodates Topcon's CV-5000 controller KB-50S as well as the CV-5000 power supply for easy printout access.

Features

- | Small footprint
- | Economical alternative refraction lane
- | Column guided cable management
- | Variety of accessories
- | Up-light room illumination
- | Right and left hand use
- | Wheelchair accessible

FS-1 accessories

OC-6



With foldable armrests



Colors

Optional accessories



Footrest

Up-down footswitch

OC-6 baseplate available with or without wheels

OC-10 / OC-12



OC-10 without armrests
OC-12 with armrests



Colors

Optional accessories



Footrest for OC-10 / OC-12

Up-down footswitch

OC-10 / OC-12 baseplate available with or without wheels


Phoropter supports



VT-670 mechanical lock arm




VT-671 balance arm

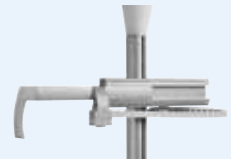


VT-72 linear arm

FS-1 optional accessories

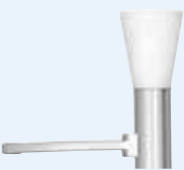

Lamp and projector bracket

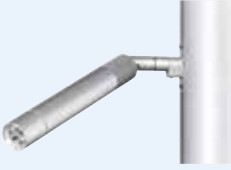

LED reading light




FS-1 dark grey cabinet

Specifications

FS-1			
Baseplate color	Grey		
Column color	Silver grey		
Cabinet color	Dark grey		
Lamp	220V, white plastic cover		
Dimensions	460 mm (W) x 500 mm (D) x 1800 mm (H)		
Weight (column & base)	53.5 kg		
Phorofter support	VT-670	VT-671	VT-72
Load capacity range	4.5 to 10 kg (adjustable)		
Stroke of phorofter arm	+ 43 inch / - 15 inch		
Max. horizontal rotation	Free		
Positioning lock	Mechanical	Counterbalanced	Autostop
Net weight of phorofter arm	6.6 kg	6.6 kg	9 kg
Linear movement	Manual		
Cable management kit	Optional		
OC Ophthalmic chairs	OC-6	OC-10	OC-12
Reclining	No		
Max seat rotation	118.5 inch	120 inch	
Foldable armrest	Yes		No
Stroke (measured at seat plane)	Min. 460 mm / Max. 560 mm		Min. 465 mm / Max. 636 mm
Footrest	Optional		
Max load capacity	175 kg		
Net weight OC	30 kg	39 kg	37 kg
Power	220V - 240V AC / 50-60 Hz		230V AC / 50-60 Hz
Electrical classification	Class I, Type B		

Dimensions

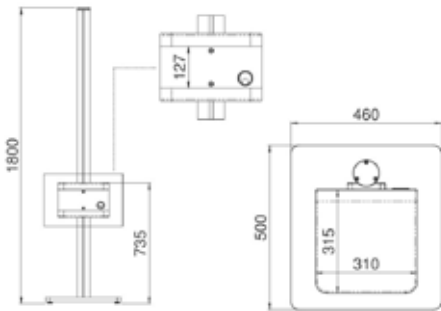

IMPORTANT

Subject to change in design and/or specifications without advanced notice.
In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.
Medical device Class I. Manufacturer: Antoni Carles S.A.



ANTONI CARLES S.A.

Topcon Europe Medical B.V.
Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O. Box 145;
2900 AC Capelle a/d IJssel; The Netherlands
Phone: +31-(0)10-4585077; Fax: +31-(0)10-4585045
E-mail: medical@topcon.eu; www.topcon-medical.eu

Topcon Danmark
Præstemarksvej 25; 4000 Roskilde, Denmark
Phone: +45-46-327500; Fax: +45-46-327555
E-mail: info@topcon.dk
www.topcon.dk

Topcon Scandinavia A.B.
Neogatan 2; P.O. Box 25; 43151 Mölndal, Sweden
Phone: +46-(0)31-7109200; Fax: +46-(0)31-7109249
E-mail: medical@topcon.se; www.topcon.se

Topcon España S.A.
HEAD OFFICE; Frederic Mompou, 4;
08960 Sant Just Desvern; Barcelona, Spain
Phone: +34-93-4734057; Fax: +34-93-4733932
E-mail: medica@topcon.es; www.topcon.es

Topcon Italy
Viale dell'Industria 60;
20037 Paderno Dugnano, (MI) Italy
Phone: +39-02-9186671; Fax: +39-02-91081091
E-mail: info@topcon.it; www.topcon.it

Topcon France Medical S.A.S.
BAT A1; 3 Route de la Révolte, 93206 Saint Denis Cedex
Phone: +33-(0)1-49212323; Fax: +33-(0)1-49212324
E-mail: topconfrance@topcon.com;
www.topcon-medical.fr

Topcon Deutschland GmbH
Hanns-Martin-Schleyer Strasse 41;
D-47877 Willich, Germany
Phone: (+49) 2154-885-0; Fax: (+49) 2154-885-177
E-mail: info@topcon-medical.de; www.topcon-medical.de

Topcon Polska Sp. z o.o.
ul. Warszawska 23; 42-470 Siewierz; Poland
Phone: +48-(0)32-670-50-45; Fax: +48-(0)32-671-34-05
www.topcon-polska.pl

Topcon Great Britain Medical Ltd.
Topcon House; Kennet Side; Bone Lane; Newbury
Berkshire RG14 5PX; United Kingdom
Phone: +44-(0)1635-551120; Fax: +44-(0)1635-551170
E-mail: medical@topcon.co.uk; www.topcon.co.uk

Topcon Ireland
Unit 276, Blanchardstown; Corporate Park 2
Ballycoolin; Dublin 15, Ireland
Phone: +353-18975900; Fax: +353-18293915
E-mail: medical@topcon.ie; www.topcon.ie



TOPCON EUROPE MEDICAL B.V.

Essebaan 11, 2908 LJ Capelle a/d IJssel, The Netherlands
Phone: +31-(0)10-4585077; Fax: +31-(0)10-4585045
E-mail: medical@topcon.eu; www.topcon-medical.eu